

BRIOSS

Cannot be
washed with
water

CHARM LIFT

ブリオス チャームリフト

BR-CA01A

BR-CA01P

BR-CA01C

PRODUCT MANUAL

Thank you for purchasing our products this time.

Please be sure to read this handling instructions carefully before using it and use it correctly.

After reading it, please take good care of this manual.



Warning

lithium-ion polymer battery built-in

- Never disassemble and take out the battery.
- Please do not disassemble it to throw it away in a dedicated collection box.

Warning

The heat-resistant temperature of lithium-ion polymer batteries is 60°C. Please do not leave it in a place where it will be hot, such as in the sun or in the car in summer.

Safety Precautions

It must be sure to protected in order to prevent damage to other people and their property.

The meaning of the general marks

The following marks show the extent of harm and damage caused by ignoring the marks and using it incorrectly.



Warning

This mark means "assumed that there is a possibility of death or serious injury."



Prohibition



Notice

This mark means "assumed that only the possibility of injury or property damage may occur".



Mandatory requirements



Warning



This product has a built-in rechargeable battery. Do not put it in the fire or heat it. Do not charge/use/leave it at high or low temperature.
It causes heat, fire, and rupture.

About the symptoms



Please do not use if you are the following people.

- People with allergies including metal allergies
- People who using in-body (planted) and wearing medical electronic devices
 - In-body implantable medical electronic equipment
 - Medical electronic equipment for life support such as artificial heart and lung
- Wearing medical electronic devices such as heart scales
- People with pre-existing conditions on the skin, such as contact dermatitis
- People who need rest or those who are not feeling well
- People with glaucoma, cataracts, and people with inflammation around the eyeball and eyes
- People with a body temperature of 38°C or higher (heat period)
example: acute inflammatory symptoms [get tired, chills, blood pressure fluctuations, etc.] are strong, and feel weak
- People whose skin is inflamed due to excessive sunburn
- People with skin problems
- People with severe spots and wounds on the skin
- People with capillary dilation
- People who can't express their intentions or people with physical disabilities
- Infants and children
- People with facial neuropathy and skin perception disorders



Please do not use on parts other than designated body parts, and the following body parts.

- Near the chest and heart
- Joints such as mucous membranes, elbows, and knees
- On top of scars and wounds after surgery, areas that are inflamed by acne and pimples
- The area where the skin is inflamed
- The part where plastic surgery was performed
- Parts where metal, plastic, silicon, etc. are embedded
- The part with a tattoo



If you are the following people, please be sure to consult with your doctor before using it.

- People with abnormal blood pressure
- People with malignant tumors
- People with heart disorders
- People immediately after giving birth from the anxiety period in the early stages of pregnancy
- People with high-level () perceptual disorders from circulatory disorders due to diabetes, etc.
- People with atopic dermatitis
- People with sensitive skin and weak skin
- People with blood circulation disorders
- People who feel abnormalities in their cranial nerves (epileptic, etc.)
- People who are being treated at medical institutions
- People who are under dental treatment (implants, orthodontics, silver teeth, etc.), people who have fillings in their teeth

About use



Do not let infants and children touch or use it while it is energized.

It will cause an unexpected accident.

Please do not apply metal or objects to the energized part.

It will cause short circuit, heat, and breakdown.

Do not use it at the same time as other beauty devices or exercise equipment.

There is a risk of harming your health.

Do not soak the body of product in water or wash it with water.

It can cause electric shock, short circuit, ignition, etc.

Do not use it continuously for a long time.

There is a risk of low temperature burns.

About "low temperature burns"

Even at a relatively low temperature (40°C to 60°C), if you touch the same place of the skin for a long time, you will get burned even if there are no subjective symptoms such as feeling hot or painful. Don't keep applying it to the same place for a long time.



Warning

About use (continued)

-  **If you feel a severe heat stimulus that seems to burn, please remove it immediately.**
Even when using it in a cold state such as in winter, it may be more sensitive to heat perception.
-  **Please make sure that the product is not damaged before use.**
Do not use if it is damaged.
-  **Please stop using it if you feel an abnormality in your body or skin, and consult a doctor in order not to harm your health and worsen your symptoms.**
It may harm your health and worsen your symptoms.
-  **If you feel strong pain or redness on your skin, please stop using it immediately.**
There is a risk of burns and injuries.

About the power supply

-  **Do not damage, process, force, pull, or twist the USB cable. Also, please don't put heavy things on or pinch them.**
It will cause electric shock, short circuit, and fire.
Do not wrap the USB cable around the main body when storing it.
The cord is overloaded and disconnected, causing a short fire, etc.
Please do not let metal or garbage adhere to the USB terminal.
It will cause electric shock, short circuit, and fire.
When cleaning, be sure to unplug the USB cable from the main unit or power supply. Also, please do not insert or remove it with wet hands.
It will cause electric shock and injury.
-  **When unplugging the USB cable, please hold the root part and pull it out without holding the cord.**
It will cause electric shock, short circuit, and ignition.

About abnormality (failure)

-  **If you find a failure, please stop using it immediately and turn off the power.**
It will cause accidents and injuries.

About use

 **Please do not use it in humid places such as the bathroom. (This product is not waterproof)**

It will cause electric shock and breakdown .

Please do not put heavy objects on the main body or let them fall.

It will cause a breakdown or an accident.

Please do not use it in the futon or with a blanket over it.

Causes of inability to breathe or suffocation

Please do not use it within 30 minutes after eating.

It will cause poor physical condition.

Do not use or store near high-temperature objects or fire.

Easy to cause deformation, malfunction or ignition.

Please do not use it for anything other than the human body.

It will cause accidents and breakdowns.

Please do not use it while wearing contact lenses.

There is a risk of eye injury.

If you use it in a state that is not moisturized enough, it may hurt your skin due to rubbing.

It will cause poor physical condition

 **When using it for the first time, please use it from mode 1 and adjust it gradually as you like.**

If you strengthen the stimulation from the beginning, it may damage your health.

Please use it within 5 minutes each time.

Continuous use for a long time will stimulate more than necessary and cause counterproductive effects and injury.

Please use it within 10 minutes in total each day.

Continuous use for a long time will cause more stings than necessary and cause counterproductive effects and injury.

Please avoid touching jewelry during use.

It will cause accidents and breakdowns.

Even if you do not use it for a long time, please charge it once every 3 months.


If you do not use it for a long time, the function will decrease and cause the failure of the rechargeable battery in the main body.

When using it, please make sure that the entire part of the energized part touches the skin.

If the contact area with the skin becomes smaller, the sting will be stronger.


Name of each part

Main body

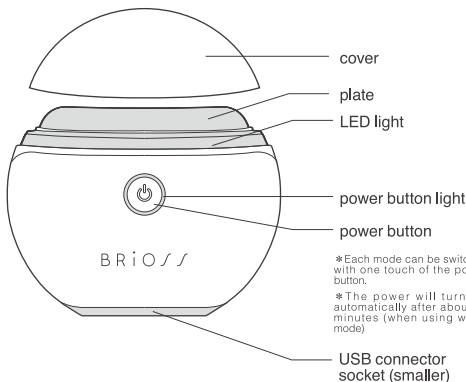

Accessories

USB charging cable

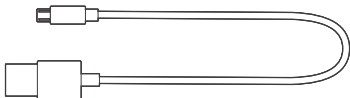


How to use

When you use it for the first time or if you don't use it for a long time, the battery may discharge and not work. Please use it after charging.

How to charge

- Please be sure to charge it with the power turned off.
- The new one is not fully charged, so please charge it before using it.

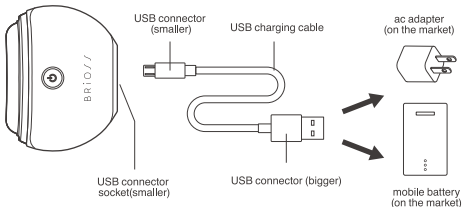

[Remark] It cannot be used while charging.

1 Insert the USB connector into each connection port.

*The power button flashes while charging.

2 Charging will be completed in about 2 hours, and the power button will turn off.

*The charging time may change depending on the condition of the battery.

*Disconnect the USB connector from each connection port.




Notice

- The recommended charging temperature is 15°C to 35°C. Outside the recommended temperature, the charging time may be longer or may not be charged.
- When you use it for the first time or when it has not been used for more than half a year, the charging time may change.
- When there is noise on the radio during charging or using, please use it in a different location.

How to use

preparation before use

Please use it after fully moisturizing with your cosmetics after washing your face.

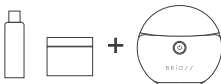

Notice

- Before use, please check the plate surface for scratches, deformation, or stains.
- This product is not waterproof. Please do not pour water on the main body.
- When using it for the first time, please use it from mode 1 and adjust it gradually according to your preference.
- Please use within 10 minutes in total each day.

* Please press and hold the power button. Each mode can be switched with one touch of the power button. After use, press and hold the power button to turn off the power.

Mode 1: Warm feeling & red LED

The red LED lights up and the plate part warms up to about 42°C. It softens your skin while gently warming it up.

Mode 2: Warm feeling & vibration & purple LED

Vibration is added to the warm mode, and the purple LED lights up. Please use it to pull up along the face line using vibration.

Mode 3: EMS & Vibration & Blue LED

* EMS mode is a micro flow, so there are individual differences in how each one feel.

EMS & vibration, blue LED lights up. The moisturizing ingredients of skin care products penetrate into the keratin liquor and provide moisturizing care.

Mode 4: Warm feeling & EMS & vibration & red-blue LED

EMS & vibration, red and blue LEDs turn on alternately. Vibration plus warmth comfortably squeezes and arranges your skin.

How to use

⚠ Notice

- If you feel strong pain or redness on your skin, please stop using it immediately. There is a risk of burns and injuries.
- If you are allergic to metal, please refrain from using it.
- EMS and the intensity of vibration may feel different, but it is not a failure.
- Please do not look directly into the LED light in LED mode.

① face line

Use **A** to slide in the direction of the arrow




② around neck and clavicle

Use **B** to slide in the direction of the arrow



③ around the forehead and cheeks

Use **A** to slide in the direction of the arrow

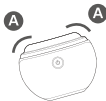

④ around eyes

Use **A** to slide in the direction of the arrow




*Do not use on eyelids

⑤ upper arm, calf

Use **B** to slide from the elbow to the armpit.



Use **B** to slide from the ankle to the back of the knee.

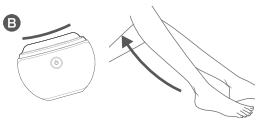

How to take care of it



警告

- Please be sure to turn off the power.
- Please be sure to remove the USB cable when cleaning the main body.
- Please wipe the main body with a damp cloth or tissue without hand washing.
- Please charge it once every 3 months.
- Please do not keep it within reach of infants and children.
- Please store away from direct sunlight, high temperature and humidity, and near moisture.
- Please untwist and store the cords.
- This product has a built-in lithium-ion rechargeable battery.
- Please be sure to follow the instructions of each prefecture and municipality for disposal.



注意

- If used with remaining moisture, it will cause injury, breakdown, and accidents.
- When disposing of the product, if there is any charge left, use the battery several times until it no longer turns on.

Is it a malfunction?

Please check according to the table below before requesting repair or service. Never disassemble or repair and remodel. There is a risk of injury or accident due to abnormal behavior.

At a time like this	How to deal with it
The power does not turn on even if you press the power button	Is it charged? (Refer to page 6)
	Did you press and hold the power button? (Refer to page 7)
The plate doesn't warm up.	Have you selected the warmth mode? (Refer to page 7)
	Isn't it out of charge? (Refer to page 6)
The power goes off during use.	Isn't it out of charge? (Refer to page 6)
	The power will automatically turn off in about 10 minutes after turning it on.(Refer to page 5)
Even if you fill it, you can only use it for a few minutes.	Have you had it more than 3 years since you purchased it?
	Please contact our customer support center.
I don't feel EMS.	EMS is a microcurrent, so there are individual differences in how you feel it.
The skin turns red or gets itchy after use.	Have you been using it for a long time?
	→Please use it while looking at the condition of your skin and body. Do not use until the redness or itching subsides.

specification

Product Name	BRiOSS Charm Lift
Model Number	BR-CA01A(blue),BR-CA01P(pink), BR-CA01C(ivory)
Power Mode	Charging
Charging Time(about)	2 hours(can be used for about 8 hours)
Rechargeable Battery	Lithium-ion polymer battery 3.7V 400mAh
Auto Off(about)	10 minutes
Power	5W
Size(about)	95×95×21mm
Weight(about)	78g
Accessories	USB charging cable, product manual
Material	Body: ABS; Head: stainless; Stopper: silicone