

How to use

preparation before use

Please use it after skin care with your cosmetics.


・We recommend to use it with cosmetics that make it easier to slide the main body so as not to put a burden on your skin.

△ Notice

- Before use, please check main body and head surface for scratches, deformation, or dirt..
- This product is not waterproof. Please do not pour water on the main body.
- Please be sure to remove accessories (necklaces, earrings, etc.).
- Please be sure to remove your makeup before using it.

turn on and select the mode and gear level.

1. Please press and hold "power/gear level changing button". Please select your favorite mode with "power/gear level changing button".
2. Press the "power/gear level changing button" and select gear level 1 to 5.
*changing gear level by short press.
*For each mode and gear level, please refer to the liquid crystal panel.

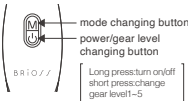

△ Notice

- When using it for the first time, please use it from gear level 1 and adjust it gradually as you like.

About vibration・EMS sensor

Turn on the power and hold the sensor(silver part) on the back of the handle, and apply the head part on the skin the EMS function will automatically activate(EMS only works in EMS mode). EMS automatically pauses when released.



HOT (warm feeling) & red LED & vibration mode

The red LED lights up, and the head part warms up to 38°C to 46°C depending on the gear level.

Green LED & vibration mode*

Please use it to pull up for each part using vibration.

EMS & Blue LED & Vibration Mode

*EMS is a microcurrent, so there are individual differences in how you feel it.

The moisturizing ingredients of skin care products penetrate into the stratum corneum of the skin to provide moisturizing care.

7

How to use

△ Notice

- If you feel strong pain or redness on your skin, please stop using it immediately. There is a risk of burns and injuries.
- If you are allergic to metal, please refrain from using it.
- EMS and the intensity of vibration may feel different, but it is not a failure.
- Please do not look directly into the LED light.

1 neck

Move it from the bottom of the neck to the chin.



2 face line

Pulls it up from chin to cheeks.



3 cheeks

Lift it from the cheeks to the bottom of the eyes.




4 around eyes

Use A to slide in the direction of the arrow




5 upper arm, calf

Move it from the elbow to the side.

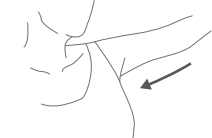

Move it from the ankle to the back of the knee.

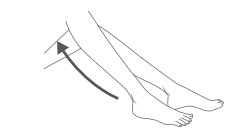

8

How to take care of it

△ Warning

- Please be sure to turn off the power.

- Please be sure to remove the USB cable when cleaning the main body.
- Please wipe the main body with a damp cloth or tissue without hand washing.
- Please charge it once every 3 months.
- Please do not keep it within reach of infants and children.
- Please store away from direct sunlight, high temperature and humidity, and near moisture.
- Please untwist and store the cords.
- This product has a built-in lithium-ion rechargeable battery.
- Please be sure to follow the instructions of each prefecture and municipality for disposal.

△ Notice

- If used with remaining moisture, it will cause injury, breakdown, and accidents.
- When disposing of the product, if there is any charge left, use the battery several times until it no longer turns on.

Is it a malfunction?

Please check according to the table below before requesting repair or service. Never disassemble or repair and remodel. There is a risk of injury or accident due to abnormal behavior.

At a time like this	How to deal with it
The power does not turn on even if you press the power button	Is it charged? (Refer to page 6) Did you press and hold the power button? (Refer to page 7)
The head doesn't warm up.	Have you selected the warmth mode? (Refer to page 7) Isn't it out of charge? (Refer to page 6)
The power goes off during use.	Isn't it out of charge? (Refer to page 6) The power will automatically turn off in about 10 minutes after turning it on.(Refer to page 5)
Even if you fill it, you can only use it for a few minutes.	Have you had it more than 3 years since you purchased it? Please contact our customer support center.
I don't feel EMS.	EMS is a microcurrent, so there are individual differences in how you feel it.
The skin turns red or gets itchy after use.	Have you been using it for a long time? →Please use it while looking at the condition of your skin and body. Do not use until the redness or itching subsides.

9

specification

Product Name	BRIOSS Neck Care Lift
Model Number	BR-NC01P/BR-NC01H
Power Mode	Charging
Charging Time(about)	2 hours
Rechargeable Battery	Lithium-ion battery 3.7V 320mAh
Auto Off(about)	5W
Power	10 minutes
Size(about)	137×78.5×37mm
Weight(about)	65g
Accessories	USB charging cable, product manual
Material	Body: ABS, Liquid Crystal: PC; Stopper: silicone

10


11

12

13

How to charge

- Please be sure to charge it with the power turned off.
- The new one is not fully charged, so please charge it before using it.

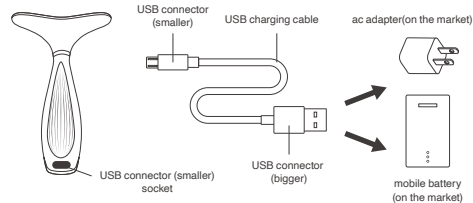

[Remark] It cannot be used while charging.

1 Insert the USB connector into each connection port.

*The power button flashes red while charging.

2 Charging will be completed in about 2 hours, and the power button light will turn off.

*The charging time may change depending on the condition of the battery.

*Disconnect the USB connector from each connection port.


△ Notice

- The recommended charging temperature is 15°C to 35°C. Outside the recommended temperature, the charging time may be longer or may not be charged.
- When you use it for the first time or when it has not been used for more than half a year, the charging time may change.
- When there is noise on the radio during charging or using, please use it in a different location.

6

Name of each part

Main body

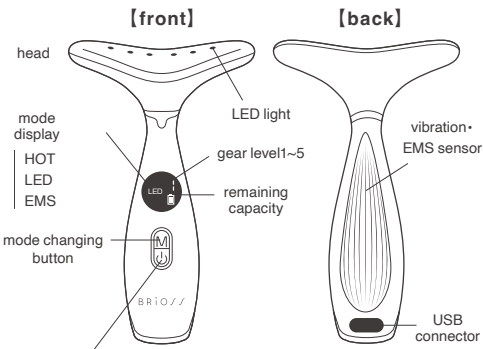


power/gear level changing button

Long press : turn on/off
short press : change gear level 1~5

- *Each mode can be switched with one touch of the power button(blue light flashes).
- *The power will turn off automatically after about 10 minutes .

Accessories

USB charging cable

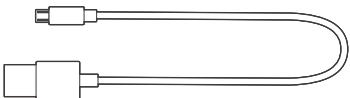

5

△ Notice

About use

- ❗ Please do not use it in humid places such as the bathroom. (This product is not waterproof)
It will cause electric shock and breakdown .
Please do not put heavy objects on the main body or let them fall.
It will cause a breakdown or an accident.
Please do not use it in the futon or with a blanket over it.
Causes of inability to breathe or suffocation
Please do not use it within 30 minutes after eating.
It will cause poor physical condition.
Do not use or store near high-temperature objects or fire.
Easy to cause deformation, malfunction or ignition.
Please do not use it for anything other than the human body.
It will cause accidents and breakdowns.
Please do not use it while wearing contact lenses.
There is a risk of eye injury.
If you use it in a state that is not moisturized enough, it may hurt your skin due to rubbing.
It will cause poor physical condition

- ❗ When using it for the first time, please use it from gear level 1 and adjust it gradually as you like.
If you strengthen the stimulation from the beginning, it may damage your health.
Please use it within 5 minutes each time.
Continuous use for a long time will stimulate more than necessary and cause counterproductive effects and injury.
Please use it within 10 minutes in total each day.
Continuous use for a long time will cause more stings than necessary and cause counterproductive effects and injury.
Please avoid touching jewelry during use.
It will cause accidents and breakdowns.
Even if you do not use it for a long time, please charge it once every 3 months.
If you do not use it for a long time, the function will decrease and cause the failure of the rechargeable battery in the main body.
When using it, please make sure that the entire part of the energized part touches the skin.
If the contact area with the skin becomes smaller, the sting will be stronger.

4

△ Warning

About use (continued)

- ❗ If you feel a severe heat stimulus that seems to burn, please remove it immediately.
Even when using it in a cold state such as in winter, it may be more sensitive to heat perception.
- ❗ Please make sure that the product is not damaged before use.
Do not use if it is damaged.
- ❗ Please stop using it if you feel an abnormality in your body or skin, and consult a doctor in order not to harm your health and worsen your symptoms.
- ❗ If you feel strong pain or redness on your skin, please stop using it immediately.
There is a risk of burns and injuries.

About the power supply

- ❗ Do not damage, process, force, pull, or twist the USB cable. Also, please don't put heavy things on or pinch them.
It will cause electric shock , short circuit, and fire.
Do not wrap the USB cable around the main body when storing it.
The cord is overloaded and disconnected, causing a short fire, etc.
Please do not let metal or garbage adhere to the USB terminal.
It will cause electric shock , short circuit, and fire.
When cleaning, be sure to unplug the USB cable from the main unit or power supply. Also, please do not insert or remove it with wet hands.
It will cause electric shock and injury.

- ❗ When unplugging the USB cable, please hold the root part and pull it out without holding the cord.
It will cause electric shock, short circuit, and ignition.

About abnormality (failure)

- ❗ If you find a failure, please stop using it immediately and turn off the power.
It will cause accidents and injuries.

3



Please do not use on parts other than designated body parts, and the following body parts.

- ・Near the chest and heart
- ・Joints such as mucous membranes, elbows, and knees
- ・On top of scars and wounds after surgery, areas that are inflamed by acne and pimples
- ・The area where the skin is inflamed
- ・The part where plastic surgery was performed
- ・Parts where metal, plastic, silicon, etc. are embedded
- ・The part with a tattoo
- ❗ If you are the following people, please be sure to consult your doctor before using it.
 - ・People with abnormal blood pressure
 - ・People with malignant tumors
 - ・People with heart disorders
 - ・People immediately after giving birth from the anxiety period in the early stages of pregnancy
 - ・People with high-level () perceptual disorders from circulatory disorders due to diabetes, etc.
 - ・People with atopic dermatitis
 - ・People with sensitive skin and weak skin
 - ・People with blood circulation disorders
 - ・People who feel abnormalities in their cranial nerves (epileptic, etc.)
 - ・People who are being treated at medical institutions
 - ・People who are under dental treatment (implants, orthodontics, silver teeth, etc.), people who have fillings in their teeth

About use



- Do not let infants and children touch or use it while it is energized.
It will cause an unexpected accident.
Please do not apply metal or objects to the energized part.
It can cause electric shock, heat, and breakdown.
Do not use it at the same time as other beauty devices or exercise equipment.
There is a risk of harming your health.
Do not soak the body of product in water or wash it with water.
It can cause electric shock, short circuit, ignition, etc.
Do not use it continuously for a long time.
There is a risk of low temperature burns.

About "low temperature burns"

Even at a relatively low temperature (40°C to 60°C), if you touch the same place of the skin for a long time, you will get burned even if there are no subjective symptoms such as feeling hot or painful. Don't keep applying it to the same place for a long time.

2

Safety Precautions

It must be sure to protected in order to prevent damage to other people and their property.

The meaning of the general marks	The following marks show the extent of harm and damage caused by ignoring the marks and using it incorrectly.
△ Warning	This mark means "assumed that there is a possibility of death or serious injury."
△ Notice	This mark means "assumed that only the possibility of injury or property damage may occur".
⊘ Prohibition	
ⓘ Mandatory requirements	

△ Warning

- ❗ This product has a built-in rechargeable battery. Do not put it in the fire or heat if. Do not charge/use/leave it at high or low temperature.
It causes heat, fire, and rupture.

About the symptoms

- ❗ Please do not use if you are the following people.
 - ・People with allergies including metal allergies
 - ・People who using in-body (planted) and wearing medical electronic devices
 - ・In-body implantable medical electronic equipment
 - ・Medical electronic equipment for life support such as artificial heart and lung
 - ・Wearing medical electronic devices such as heart scales
 - ・People with pre-existing conditions on the skin, such as contact dermatitis
 - ・People who need rest or those who are not feeling well
 - ・People with glaucoma, cataracts, and people with inflammation around the eyeball and eyes
 - ・People with a body temperature of 38°C or higher (heat period)
 - ・example: acute inflammatory symptoms (get tired, chills, blood pressure fluctuations, etc.) are strong, and feel weak
 - ・People whose skin is inflamed due to excessive sunburn
 - ・People with skin problems
 - ・People with severe spots and wounds on the skin
 - ・People with capillary dilation
 - ・People who can't express their intentions or people with physical disabilities
 - ・Infants and children
 - ・People with facial neuropathy and skin perception disorders

1

BRIOSS

NECK CARE LIFT

ブリオス ネックケアリフト

BR-NC01P

BR-NC01H


Cannot be washed with water

PRODUCT MANUAL

Thank you for purchasing our products this time.

Please be sure to read this handling instructions carefully before using it and use it correctly.

After reading it, please take good care of this manual.

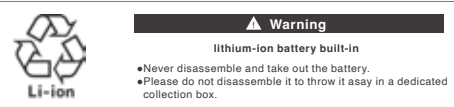

△ Warning

The heat-resistant temperature of lithium-ion polymer batteries is 60°C. Please do not leave it in a place where it will be hot, such as in the sun or in the car in summer.